

PRODUCT INFORMATION

Common Name SGN-70, h1F6

Synonyms CD70;CD27L;CD27LG;TNFSF7

Applications ELISA; Flow Cyt

Recommended

Dilutions

ELISA 1:5000-10000; Flow Cyt 1:100

Lyophilized from sterile PBS, pH 7.4. Normally 5 % – 8% trehalose is added as protectants before lyophilization. Please see Certificate of Analysis Formulation & Reconstitution

for specific instructions of reconstitution.

Host Species Homo sapiens

IgG type lgG1 Reactivity Human **CD70 Target Uniprot ID** P32970

Description Anti-CD70 (vorsetuzumab biosimilar) mAb

Delivery In Stock

> Store at -20°C to -80°C for 12 months in lyophilized form. After reconstitution, if not

intended for use within a month, aliquot and store

Storage & Shipping at -80°C (Avoid repeated freezing and thawing). Lyophilized proteins are shipped at ambient

temperature.

Research grade biosimilar. Not for use in

Background therapeutic or diagnostic procedures for humans

or animals.

Usage Research use only

Conjugate Unconjugated

All DIMA recombinant antibodies are genuinely

generated by DIMA Biotech. They are all under patent application. Any protein sequencing or reverse engineering attempt is prohibited. We are **DIMA Disclaimer**

actively scrutinizing all patent application to ensure no IP infringement.

Email: info@dimabio.com Website: www.dimabio.com





Anti-CD70 (vorsetuzumab biosimilar) mAb ELISA

0.2µg of Human CD70, mFc-His Tagged protein per well

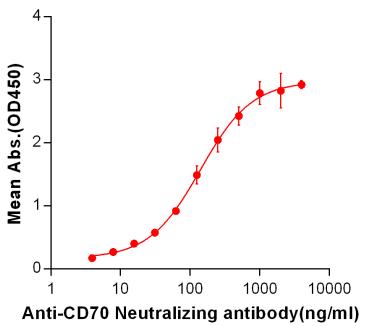


Figure 1 . ELISA plate pre-coated by 2 μ g/ml (100 μ l/well) Human CD70, mFc-His tagged protein PME100042 can bind Anti-CD70 Neutralizing antibody in a linear range of 3.91-139.9 ng/ml.

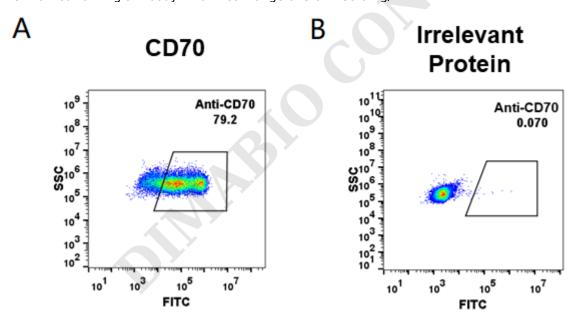


Figure 2. HEK293 cell line transfected with irrelevant protein (B) and human CD70 (A) were surface stained with anti-CD70 neutralizing antibody 1µg/ml (vorsetuzumab) followed by Alexa 488-conjugated anti-human IgG secondary antibody.

Email: info@dimabio.com Website: www.dimabio.com





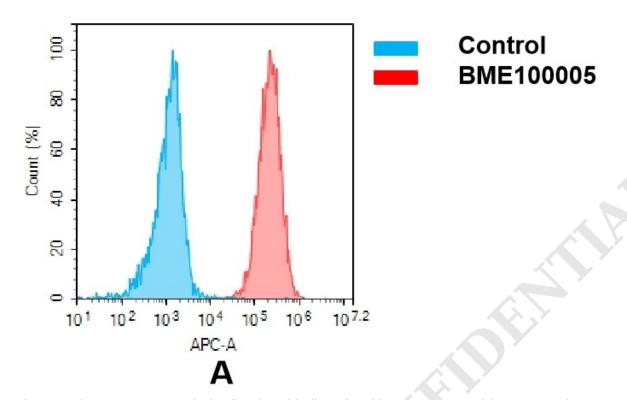


Figure 3. Flow cytometry analysis of antigen binding of anti-human CD70 mAb(BME100005). (A) A clear peak shift of BME100005 was seen compared to the control when incubated with CD70-expressing THP1 cells, indicating strong binding of BME100005 to CD70. Antibodies were incubated at 2 μ g/mL.



