

## **PRODUCT INFORMATION**

**Applications** Antibody internalization labeling reagent

**Detection method** Cell viability detection with MTT, CCK8, or CTG

The product has a MW of 34 kDa **Molecular Weight** 

Lyophilized from sterile PBS, pH 7.4. Normally 5 % Formulation & 8% trehalose is added as protectants before Reconstitution lyophilization. Please see Certificate of Analysis

for specific instructions of reconstitution.

The DiTagTM MMAE IgG labeling reagents can be used for human IgG1, IgG2 and IgG4, rabbit IgG, mouse IgG2a and IgG2b. IgG type

Recommended We recommend test antibody to mix with

**Dilutions** AME100003 at 2:1 in molar

Description DiTagTM MMAE IgG labeling reagent

**Delivery** in Stock

**Background** 

The reagents are supplied in lyophilized form. We

recommend storing the vial(s) at -20°C,

Storage & Shipping desiccated and protected from light. Once

reconstituted, the reagents can be stored at 2-8°C for 1~2 weeks, or with 50% glycerol at -20°C. DiTag™ MMAE IgG labeling reagents provide an

easy solution for quantifying antibody internalization activities. Leveraging VcMMAE (MC-VC-PABC-MMAE) conjugated to an Fc binding protein, these reagents bind to IgG antibodies from various species, resulting in the formation of a VcMMAE-labeled antibody-reagent complex.

Upon antibody internalization, the cleavable linker MC-VC-PABC is enzymatically cleaved by cathepsin B, a protein overexpressed in multiple cancer types. This enzymatic cleavage triggers the release of PABC-substituted MMAE, forming an unstable intermediate that liberates the free drug. Measurement of cell killing or inhibition allows researchers to evaluate the efficiency of antibody internalization into cells. This critical information enhances our understanding of the cellular uptake mechanism of antibodies and aids

> Email: info@dimabio.com Website: www.dimabio.com

in assessing their efficacy in targeted therapies or diagnostic applications.

Usage Research use only

Address: Wuhan institute of Biotechnology B7, Biolake No.666 Gaoxin Road, Wuhan, Hubei, China Telephone: +1 2409940618(USA) /+86-18062749453(China) /+86-400-006-0995(China)





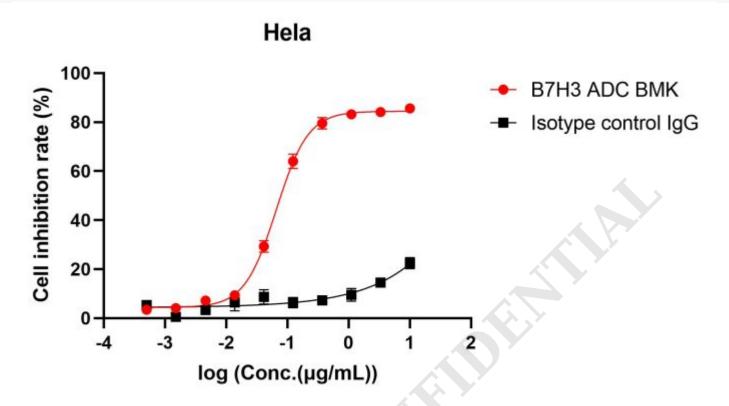


Figure 1. Cell inhibition rate of Hela detected by CCK8 method. The IC50 of B7H3 ADC BMK is 66ng/ml, indicating specific internalization.

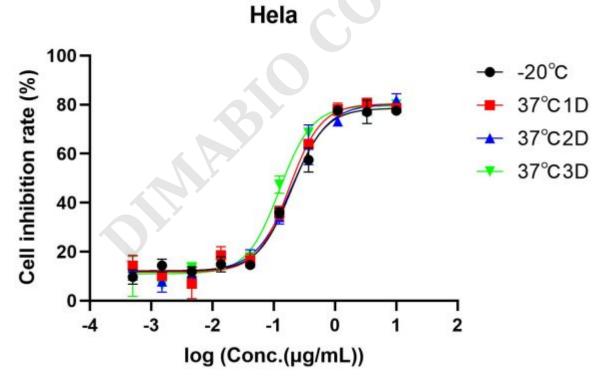


Figure 2. Accelerated stability test of AME100003. Following lyophilization, samples were stored at -20°C (black), 37°C for 1 day (red), 37°C for 2 days (blue), and 37°C for 3 days (green), separately. Upon reconstitution, the cell inhibition rate of each sample was determined using the CCK8 method. The data indicate excellent stability for all samples.

Address: Wuhan institute of Biotechnology B7, Biolake No.666 Gaoxin Road, Wuhan, Hubei, China Telephone: +1 2409940618(USA) /+86-18062749453(China) /+86-400-006-0995(China)

Email: info@dimabio.com Website: www.dimabio.com

