

## PRODUCT INFORMATION

<b>Tag</b>	<b>Warning:</b> Undefined variable \$hasAttributeValueDescription in C:\wwwroot\mirror\dimabio.com\wp-content\plugins\woocommerce-print-products\publicclass\woocommerce-print-products-public.php on line 2806 C-Flag&Strep Tag
<b>Target</b>	PACR
<b>Synonyms</b>	PAC1; PAC1R; PACAPR; PACAPRI
<b>Description</b>	Human PACR-Strep full length protein-synthetic nanodisc
<b>Delivery</b>	In Stock
<b>Uniprot ID</b>	P41586
<b>Expression Host</b>	HEK293
<b>Protein Families</b>	Druggable Genome, GPCR, Transmembrane
<b>Protein Pathways</b>	Neuroactive ligand-receptor interaction
<b>Molecular Weight</b>	The human full length PACR-Strep protein has a MW of 53.3 kDa
<b>Formulation &amp; Reconstitution</b>	Lyophilized from nanodisc solubilization buffer (20 mM Tris-HCl, 150 mM NaCl, pH 8.0). Normally 5% - 8% trehalose is added as protectants before lyophilization. Please see Certificate of Analysis for specific instructions of reconstitution.
<b>Yefei_Storage</b>	Store at -20°C to -80°C for 12 months in lyophilized form. After reconstitution, if not intended for use within a month, aliquot and store at -80°C (Avoid repeated freezing and thawing). Lyophilized proteins are shipped at ambient temperature.
<b>Background</b>	Type I adenylate cyclase activating polypeptide receptor, which is a membrane-associated protein and shares significant homology with members of the glucagon/secretin receptor family. This receptor mediates diverse biological actions of adenylate cyclase activating polypeptide 1 and is positively coupled to adenylate cyclase.
<b>Usage</b>	Research use only
<b>Conjugate</b>	Unconjugated

### ELISA assay to evaluate PACR-Strep-Nanodisc 0.2µg Human PACR-Strep-Nanodisc per well

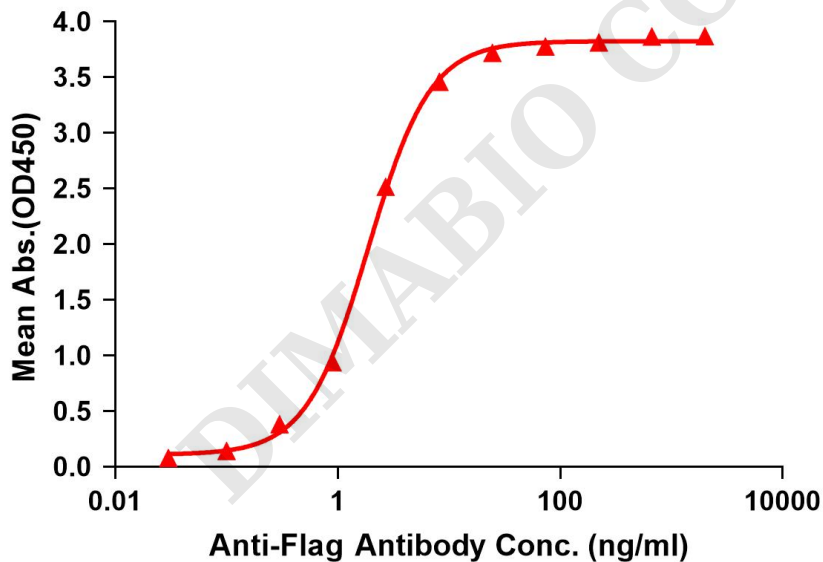


Figure 1. Elisa plates were pre-coated with C-Flag&Strep Tag PACR-Strep-Nanodisc (0.2µg/per well). Serial diluted anti-Flag monoclonal antibody solutions were added, washed, and incubated with secondary antibody before Elisa reading. From above data, the EC50 for anti-Flag monoclonal antibody binding with PACR-Strep-nanodisc is 1.887ng/ml.



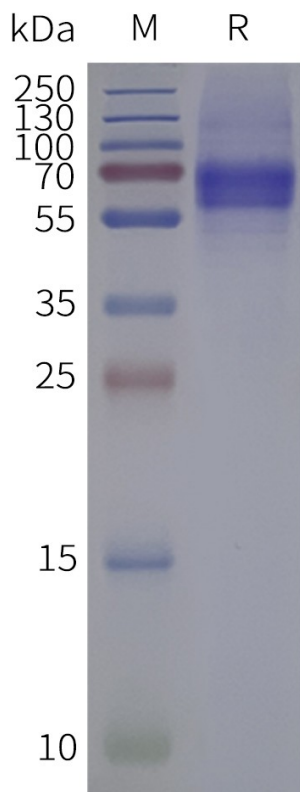


Figure 2. Human PACR-Strep-Nanodisc, C-Flag&Strep Tag on SDS-PAGE

DIMABIO CONFIDENTIAL

